


# PADDLE BLENDER HOMOGENISER



HBM-400G

This instrument specialises in the sterile homogenisation of samples. The gentle yet rapid blending action minimises cell and tissue damage whilst maintaining efficiency.

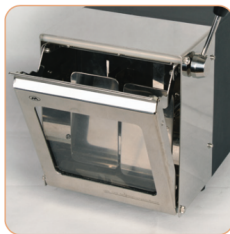

## Easy Maintenance

Front loading door is removable, and only stainless steel makes contact with the sample bag, making the whole machine easy to clean and maintain



## Touch-screen display

Adjustable blending time and frequency. Shows parameters and current conditions of the machine. Memory for 6 different settings, making experiments repeatable.



## Sterile bag

Reduces risk of sample cross-contamination.



HBM-400B

### Included:

Paddle Blender Homogeniser x1, Homogeniser Bags x25, Rod Clamp File x5, Water Catching Plate x1, Power Cable x1, Instruction Manual x1

Series	HBM	
Models	HBM-400B	HBM-400G
Control Mode	Touch-screen	
Blender Distance Adjustment	Manual	Auto (touch-screen)
Digital Timer Range	10, 30, 60, 90, 120, 180, 600s, continuous	0-60min, continuous
Speed Range	6-9 strokes/sec	4-10 strokes/sec
Memory Capacity	6 presets	
Sample Volume	50-400mL	
Leakage Alarm	-	✓
Weight (kg)	16.5	
Dimensions (mm)	410x270x290	
Sterile Bag Dimensions (mm)	300x185	